



permobil



## JUST THE FACTS

Combining stability and agility, the C300 Corpus® 3G gives you the high caliber performance you expect from Permobil. Its compact size lets you easily navigate tight spaces, while the powerful motors take you anywhere you need to go. The Corpus 3G Ergo Backrest mimics the shape and comfort of today's Corpus Backrest, but provides more functionality and adjustment. A powerful center-mounted actuator offers up to 50° of tilt, and users can recline up to a full 175°. The legrests are curved for better support, and the armrests feature built-in channels for cable routing, which significantly reduce the need for messy wire ties. All of this evidence of the ideal power wheelchair...case closed.


**C300**  
CORPUS 3G

**POWER FUNCTIONS**




Power tilt and recline with elevating legrests

**SHROUD COLORS**

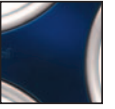


Crystal Blue

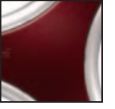


Vintage Wine



Metallic Green

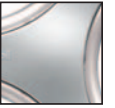

Platinum Silver

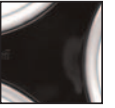


Black Onyx



Chili Pepper Red



Popstar Pink

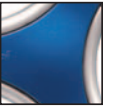

Titan Blue



50° of tilt with two range selections



175° power recline




8" power seat elevator



Power adjustable legrests



Independently adjustable legrest height

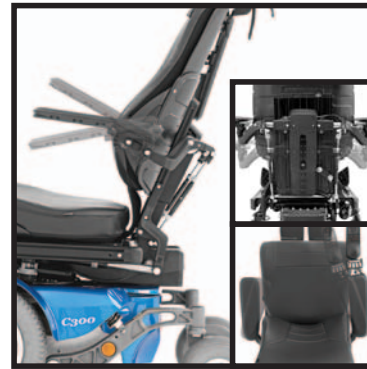

Independently adjustable armrest height, angle and pitch

**SPECIFICATIONS**

Maximum Speed	6.5 MPH (High Speed Package)	Seat Widths	17", 19", 21", 23" (Adjustable)
Driving Range	15-19 Miles	Seat Depths	14" - 22" (Adj. in 1" increments)
Max User Weight	300 lbs	Backrest Heights	20", 23" - 28"
Ground Clearance	3"	Seat To Floor Height	17.5"
Battery Capacity	Group 34	Power Tilt	50°
Base Width	24"	Power Recline	85° - 175°
Length	43"	Mechanical Shear Reduction	5.25"
Recommended Battery	Group 34	Power Legrest	90° - 170°
Gradient Capacity	10° Up 6° Dn	Legrest Articulation Amount	8"
Power Seat Elevator	8"	Electronics (standard)	R-Net
Arm Pad Options	10", 13", 16"	Alternative Drive Controls	Yes w/R-Net



THE POWER OF MOBILITY

**Permobil Inc.**  
**300 Duke Drive**  
**Lebanon, Tennessee 37090**  
**800.736.0925**  
**www.permobil.com**